

PCL Separation Science Lab

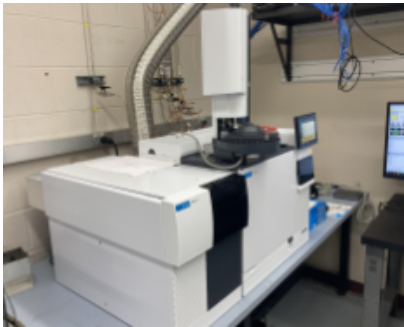
Principle

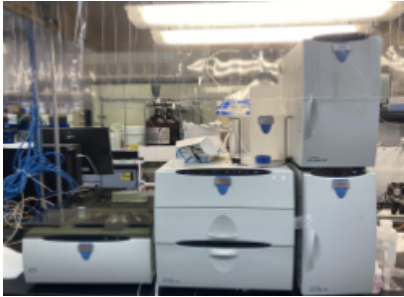
Separation Science

Capabilities

- GC, GC-MS
- HPLC
- IC
- GPC (SEC)
- UPLC-MS TOF with ASAP
- TDGC-MS
- Preparative Chromatography

Equipment

	Equipment Type	Information
	GC-MS/MS	Purity Analysis and Trace Contamination GC • 9 dual channel GCs • Flame Ionization Detector (FID) • Thermal Conductivity Detector (TCD) • Electron Capture detector (ECD) HPLC/UHPLC • 3 LC systems • UV/Vis • Charged Aerosol Detector (CAD) GC-MS/MS (QQQ) • ppb level detection
	Headspace GC/GCMS	Unknown Volatiles & Decomposition GC-MS • 3 GC-MS systems • 2 multidetector systems (FID/MSD): simultaneous collection of FID signal and MSD signal for both qualitative and quantitative analysis TD/Py GC-MS • Product deformation • Evolved Gas Analysis (EGA): 50°C - 1000°C • Identification of unknown base resins • Complimentary to TGA • Thermal Desorption/Pyrolysis (TD/Py): 50°C - 1000°C • Identification of off-gassing at various temperatures Headspace GC/GCMS • Analysis of volatiles in solid and liquid samples: Ambient - 300° • FID and MSD detection for qualitative and quantitative measurements

	Ion Chromatography	Ion Quantification IC - Quantification of residual ionic species that are extractable into aqueous medium - Chloride, fluoride, sulfide, acetate, methane sulfonate, etc. Combustion IC (CIC) - Quantification of total ionic species that are generated upon combustion of polymeric material - Typically chloride and fluoride
Multi-Detector GPC  HTGPC	HTGPC	Molecular Weight Distribution GPC 6 instruments each dedicated to a polymer type (mostly RI detection, but some UV/Vis) High Temp GPC (HTGPC) - For chemically resistant materials - Capable of temps from ambient up to 220°C
	Multi-Detector GPC	Multi-Detector GPC - Used for non-routine analyses - Dual angle LS, RI, and Viscometry - Can be used for traditional GPC and also for measurement of absolute molecular weight and intrinsic viscosity to gain insight into conformational characteristics and other properties of polymers
Mission	Separation Science mission	
Techniques	Separation Science	

Team

Separation Science Lab

	Contact	Ext	Office
	Scott Sudbury Research Scientist Separation Science	M: 404-274-4134	Office N2083 Lab N2102, N2106
	Cheston Husein Separation Science	770-772-8419	Lab N2104 Cubical N2048.3
	Brayden Koster Separation Science & Spectroscopy	770-772-8361	Lab N2104
	Steve Myers Separation Science	770-772-8696	Lab N2104
	Olga Mytykh GC, GC-MS, IC	770-772-8665	Lab N2102
	Daniel Quiroz GC, GC-MS, IC	770-772-8663	Lab N2102 Cubical N2048.4